




# YA-D5-9

ELECTRIC HOSPITAL BED

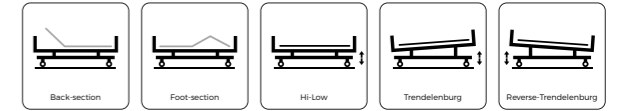

**PRODUCT  
DETAILS**



**Remote control**

The integrated remote control with clear icons allows nurses to quickly and clearly understand the sections they need to control.



**Soft link**

Soft link makes platform more smooth to let the patient feel more comfortable when adjust the position



**Motors**

CE mark motors can be replaced according to requirements. Strong stability, great support without any noise.



**Detachable Headboard**

The lock structure at the end of the bed is convenient for quick removal in an emergency, allowing medical personnel to carry out rescue.



**PP safety railing**

Tuck-away PP safety railing with damper for slow dropping, special patented design to avoid accidental falling from patient's on bed operation



**Central locking system**

Central locking system with dual sided wheels is with excellent stability and convenient for medical staff to manipulate

- Four split side rail, drop down under the bed platform
- Hand remote control
- Detachable head/foot board with easy lock
- With articulated ABS soft joint between different platform sections
- Central brake system with wheel diameter of 125mm
- Auto-counter to provide the patient from sliding onto foot end
- Manual CPR on both sides(optional)

1	Overall size(H*W*D)	2160*1050 *450-700mm
2	Lying size	1910*900 mm
3	Backrest adjustment	0-75°
4	Footrest adjustment	0-40°
5	Height adjustment	450-700mm
6	Trendeleburg	0-12°
7	Revere- Trendelenburg	0-12°
8	Safe loading capacity	205kg

