




# MC-D05

GYNECOLOGICAL TABLE

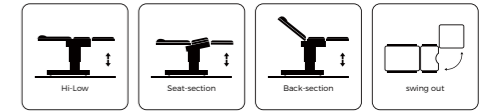

  
**PRODUCT  
DETAILS**



### Castor

The table can be moved easily or in a stable position thanks to the special leg elevation system.



### Mattress

Modular foaming high-density mattress makes people more comfortable, easy for hygiene



### Hand / Foot control

Sometimes, the doctor needs the table to be adjusted when a certain hand operation is being done, the foot pedal switch can free the hand to adjust the table which is a vast convenience.



### Foot section

The foot section can be rotated to one side and removed for easy delivery and exam



### Column and Base

The main structure, specially for column and base cover, is made of high quality stainless steel which is resistant to corrosion.



### Armrest lock

T-shaped unique armrest lock, can rotate different angles to suitable for patient

- Unique profile ergonomic design which is beautiful and convenient.
- All base and column cover is made of 304 medical stainless steel.
- Mattress adopts high density sponge and is once foam molding, with no seam.
- Hand control and foot pedal, mutually independent and is easily for operation.
- Equipped with hidden sewage basin.
- Armrest swing out  $\geq 90^\circ$

1	Model	MC-D05
2	Overall size	1780* 600 * 680-880mm
3	Back section adjustment	-7°-45°
4	Height adjustment	680-880mm
5	Seat adjustment	0°-30°
6	Loading capacity	225kg
7	swing out	0-90°

