



# OVERBED TABLE

MODEL: YA-T05

# YA-T05

YA-T05 over bed table is with brilliant features of the wood-grain laminate top plate, chrome painting frame making the bed table with an elegant appearance, spring adjusted height for convenient height operation, unique patented adjusting system can lift the table automatically up from underneath force which will avoid the accidental turning.





ITEM	PARAMETERS
Overall size	L710*W370*H700-1140mm
Table top size	710*370mm
Thickness of table top board	20mm
Castor	4 x 50mm (2 w/brake)
Safe working load	25kg

**01 TOP BOARD**  
Wood-grain laminate top with molded edge for damage protection

**02 PROTECTIVE DEVICE**  
After installation, remove the black screws to adjust the height

**03 CASTOR**  
High quality medical caster 2 inch mute wheel,2 with brake

**04 SPRING**  
Air spring for height adjustment lifting700-1140mm

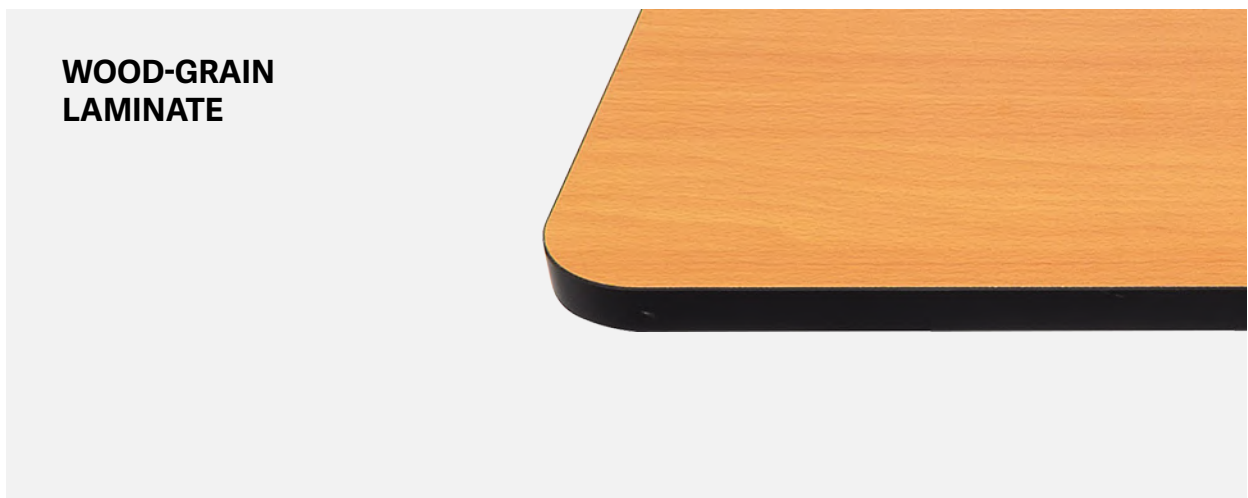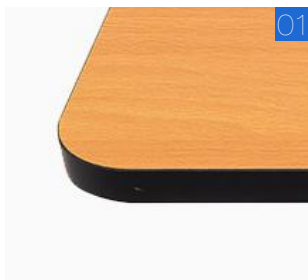



You Tube



Medik

Tel : +86 15250358182 / Fax : +86 51258688450

E-mail : [info@medik-medical.com](mailto:info@medik-medical.com)

Address : No.3 Jinxiu Road,Jinfeng Industrial Zone,Jinfeng  
Town,Zhangjiagang,Jiangsu Province,China