

YA-DH5-2U

Home Care Bed



*Better care
Better comfort!*

YA-DH5-2U

Home Care Bed

DESCRIPTION

YA-DH5-2U is an elegant nursing bed designed for home care and nursing center. Design with Integration of patient's safety, comfort, environment harmony, the nursing bed is very well recognized with the market. Warm brown wood grain design increases the patient's sense of belonging and makes the atmosphere more warm. Further, thanks to the special packing design, the bed is also welcome with private market due to minimized and optimized packages.

FEATURES

- Electrically operated profiling bed for community and nursing home use
- Wooden head and footboard with aluminum columns around
- Adjustable self-locking side rails providing additional safety which can be released at both head and foot ends
- Safe work load (SWL) 180Kg
- Variable height (400mm - 800mm)
- Remote handset control
- High-quality castors (castor diameter: 100mm)
- Quick and easy assembly in limited space

STANDARD SPECIFICATION

Model	YA-DH5-2U
Overall size	2160*1030 *400-800mm
Lying size	1900*900 mm
Backrest adjustable	0-75°
Footrest adjustable	0-40°
Height adjustable	400-800mm
Trendelenburg	0-12°
Reverse-Trendelenburg	0-12°
Loading capacity	180kg



CONTROL SYSTEM

The bed control system uses the German DEWERT, and the linear heavy duty type actuator guarantees the safe load capacity will be 180 KG weight and can adjust the different positions .



AUTO COUNTER FUNCTION

The bed can be adjusted with one button to a comfortable seating position to watch TV or read book.




BED FRAME

Four pieces of welded bed board, light but strong, with mattress holders.



SIDE RAIL

The full-length side rail minimize the risk of the patient falling down, and the height of the guardrail can be adjusted by pressing button.

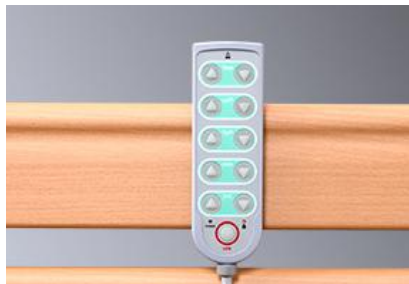

REMOTE CONTROL

The hand control with clear icon is beneficial for medical staff which can control the bed position from a remote distance.



CALF BOARD

Calf board can be adjusted by racks, which can raise level the feet, promote blood circulation if necessary.



CASTOR

Four ϕ 10 cm swivel casters with brakes , easy to move the patient and fixed in any desired position.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US