

# YA-DH5-3

Home Care Bed



*Better care  
Better comfort!*

# YA-DH5-3

# Home Care Bed

## DESCRIPTION

This YA-DH5-3 is extra-low home care bed with a minimum height of 300mm, the low height will minimize the consequence from bed falling. Full-length guardrails gives a perfect protection for patient from falling off. The entire wooden design generates the bed with an home style bed environment.

## FEATURES

- Electrically operated profiling bed for community and nursing home use
- Safe work load (SWL) 180Kg
- Variable height (400mm - 800mm)
- Wooden head and footboard with aluminum columns around
- Adjustable self-locking side rails providing additional safety which can be released at both head and foot ends
- Remote handset control
- High-quality castors (castor diameter: 100mm)
- Quick and easy assembly in limited space

## STANDARD SPECIFICATION

Model	YA-DH5-3
Overall size	2160*1030*300-700mm
Lying size	1900*900mm
Backrest adjustable	0-75°
Footrest adjustable	0-40°
Height adjustable	300-700mm
Trendelenburg	0-12°
Reverse-Trendelenburg	0-12°
Loading capacity	180kg



**GERMANY DEWERT CONTROL SYSTEM**

Total 4pcs linear heavy duty type actuator will guarantee different position adjustment and safety loading 180kg.



**BED PLATFORM**

Sturdy metal welded platform surface with mattress holder, 4 section board.



**FULL LENGTH SIDE RAIL**

The full-length side rail prevents the patient from falling down, and the height of the guardrail can be adjusted by pressing button.



**SIDE RAILING STOPPER**

There are stoppers under the headboard and footboard, which can prevent the side rail from sliding down and prevent safety accidents.



**PRECISION LIFTING COLUMN**

4pcs precision seamless lifting columns will guarantee bed adjustment noiseless.



**U-SHAPED FRAME**

The special U-shaped frame is closer to the ground, so that the minimum height of the home care bed is 300mm.




**EXTRA LOW HEIGHT**

The extra low position is designed to alleviate the falling consequence from severe injury.



**HIGH POSITION**

Max height 700mm is a comfortable height for doctor to examine without straining the doctor back.

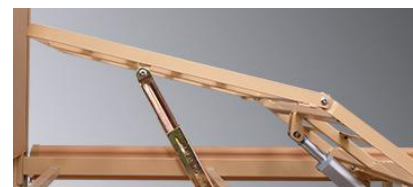

**CALF BOARD**

Calf board can be adjusted by racks, which can raise level the feet, promote blood circulation if necessary.



**REMOTE CONTROL**

Hand Remote control with electric CPR function, millions of secondary button service life, clear function identification, easy operation.



**CASTOR**

Unique four  $\Phi 100\text{mm}$  swivel casters, high support, impact resistance, not easy to break.

## ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**info@medik-medical.com | +86 139 6223 1292**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US