

YA-AS09

Non-Magnetic MRI Stretchers



YA-AS09 Non-Magnetic MRI Stretcher is suitable for use in 1.5t and 3.0t MRI examinations. In a strong magnetic environment, there will be no stretcher movement or transfer, no suction, no interference on the results of MR examinations, and convenient transfer of patients. The stretcher can be separated, and the board can be separated from the frame for quick transfer.

Technical Specification

ITEM	PARAMETERS
Size	190*55*40-88cm
Load capacity	160KG
Weight	27KG


Main features

- The frame is composed of a bed surface and a car body. The bed surface can be separated from the car body.
- The stretcher is equipped with self-locking side rails
- After strict demagnetization treatment, aluminum alloy parts will be non-magnetic
- The backrest is adjusted by the handle
- All fixing parts are demagnetized screws



Fixed slot

With four fixed slots, it can be placed quickly while effectively preventing the movement of the bed surface.

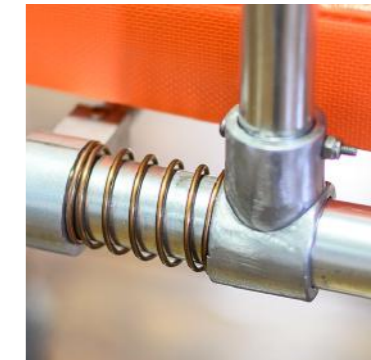

Self-locking side rail

The Self-locking side rail can be retracted quickly and is convenient for the operator to use.



Back adjust

Back is adjustable which makes the wounded feel comfortable. Back rest lifting 0-60°



IV stand

The infusion stand is also non-magnetic, and the infusion stand can be stored on one side of the stretcher. There are infusion holes on both sides of the stretcher.



Height adjust

The end of the stretcher is equipped with a handle that can adjust the height, and the "X"-shaped.



Castor

5' silence non-magnetic plastic wheel with PU ring ,advanced protection universal brake wheel,move easily , stable and durable.

