

YA-PS15

YA-PS15 is a manual patient transfer stretcher with height and back adjustment. It can safely transfer the patient from the stretcher to the bed. innovative side-rail design allows medical staff to carry out this action efficiently.



Back section

The back section is adjusted via gas spring to any wanted degree between 0~70



Height adjustment by crank

The height adjustment is controlled by crank, which is convenient for nurses to control




Oxygen cylinder holder

The oxygen cylinder holder is specially designed under the backboard of the stretcher frame



IV Pole

The storage hole of the iv pole is hidden under the back board

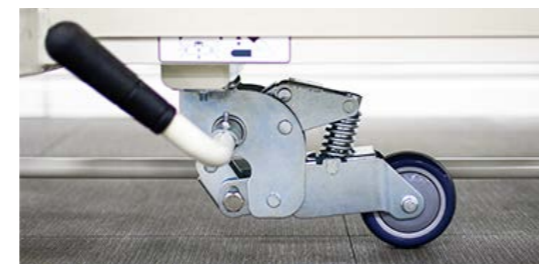

5th castor

With 5th retractable wheel in the middle for direction



Central brake castor

Diameter 200mm big castor with central control brake has good passability during transportation



BASE



URINE BAG HOOK



IV POLE HOLDER

Overall size

1930*640*520-820mm

Back section adjustment

0-70°

Height adjustment

520-820mm

Caster

Diameter 150mm

Loading Capacity

250kg

Main features

- Backrest adjustable in a pneumatic way
- Manual luxurious stretcher
- ABS tuck-away railing, collapsible
- With 5th retractable wheel in the middle for direction
- Central brake system with 150mm diameter wheel
- ABS platform easy for cleaning
- Manual height adjustable

Accessories

- Mattress
- IV pole

