

# YA-SS02

Lightweight Stair Stretcher



*Better care  
Better comfort!*

# YA-SS02

# Lightweight Stair Stretcher

## DESCRIPTION


The YA-SS02 is a light weight stair stretcher made of high strength aluminum alloy. It is widely used in the home, school, hospitals etc for those who can not move alone during emergency evacuation specially through stairs, which is very easy to store and transport with its light weight and folded structure.

## FEATURES

- This stair stretcher is made of high-strength aluminum alloy materials.
- It is an ideal facility for emergency staff to transfer patient to go up and down stair.
- Four wheels enable the patients to be moved over most surfaces with minimal effort.
- The stretcher approves the patient safety during transport process because of its 2pieces of restraint straps with quick-release buckle.
- Under the front of the stretcher, it has two flexible uplift rods.
- The structure can be folded and it saves the space in ambulance car.
- Light-weight and easy for sterilization and clean.

## STANDARD SPECIFICATION

Size (Unfolded )	900*510*910mm
Size (Folded)	900*170*590mm
Loading capacity	159kg
Carton Size (1pc/carton)	960 * 160 * 540mm
N.W.	8kg
G.W.	9.2kg

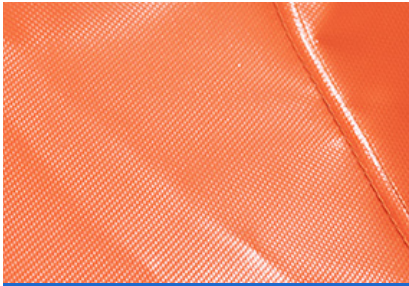

**Material**

Coated Oxford cloth is tough and not easy to damage, and the surface can be easily cleaned of blood stains and stains.



**Castors**

There are 4 castors which make it easier to transport the injured on a flat ground, and 2 with brakes, which can be fixed at designated positions.



**Foldable design**

The whole stretcher include the handle can all be folded, which can reduce the volume and put it in any desired position.



**Folding handle**

Stair stretcher can be folded to the smallest volume in one step for easy storage.



**Carrying handle**

Adjustable carrying handles can be adjusted in length and angle for maximum convenience when going up and down stairs.

## ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**info@medik-medical.com | +86 139 6223 1292**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.