

MC-C08

Gynecology table



*Better care
Better comfort!*

MC-C08

GYNECOLOGY TABLE

DESCRIPTION


MC-C08 simple gynecological exam table is made of high-quality stainless steel round tubes welded together. It is sturdy, strong and easy to clean, and can meet the needs of most gynecological examinations.

MAIN FEATURES


- Gynecological table with high quality stainless steel frame.
- Three sections with independent controls.
- Mounted on 4 adjustable rubber feet.
- Back section adjustable by lifting bar.
- Pelvic and foot section adjustable by ratchet.
- Fitted with patient hand grips.
- leg support holders adjustable in height and width.
- Fitted with removable filth basin.

TECHNICAL PARAMETER

Model	MC-C08
Overall Size	1800*800*H750mm
Back section adjust	0-60°
Pelvic section adjust	0-10°
Leg section adjust	-90°-0°
Mattress thick	50 mm
Loading capacity	180kgs

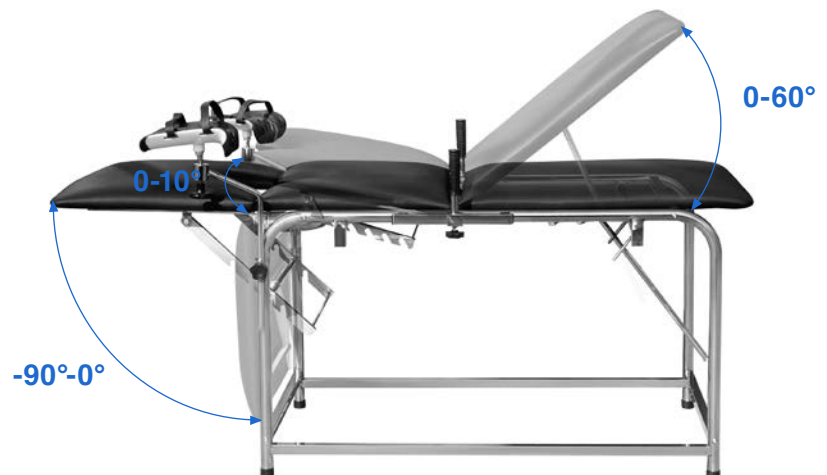

MAIN DETAILS

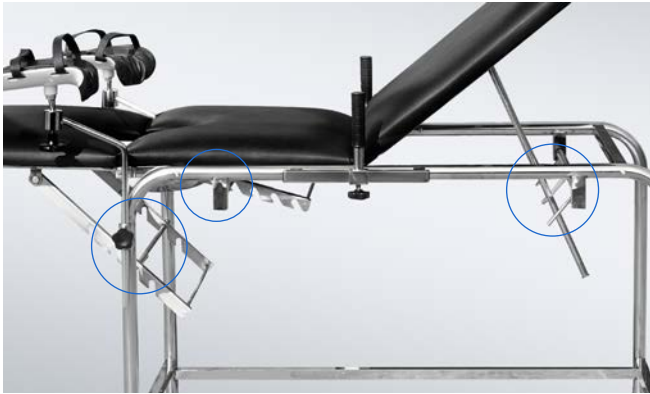

Adjustment mechanism

The back, Pelvic and leg part of the gynecological exam table are adjusted through 3 independent ratchet mechanisms.




Dirt Basin

Stainless steel dirt basin is placed under the semicircular hole of the gynecological table and can be used to collect fluids.




Rubber feet

The 4 legs of the gynecological table are equipped with non-slip rubber feet, which can be stably fixed on the ground without moving.

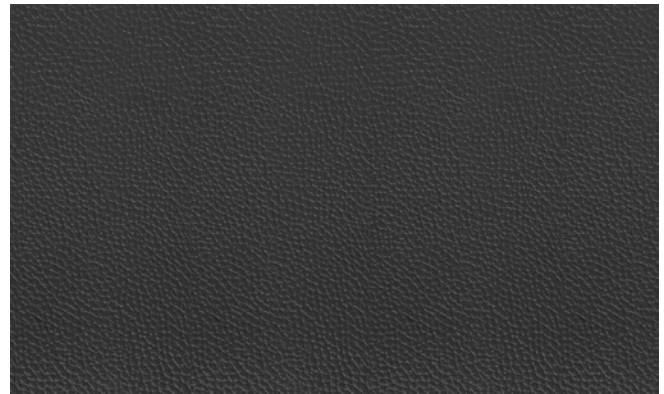

Foam and leather

High-density sponge with thickness of 5cm, cover leather is waterproof and wear-resistant.

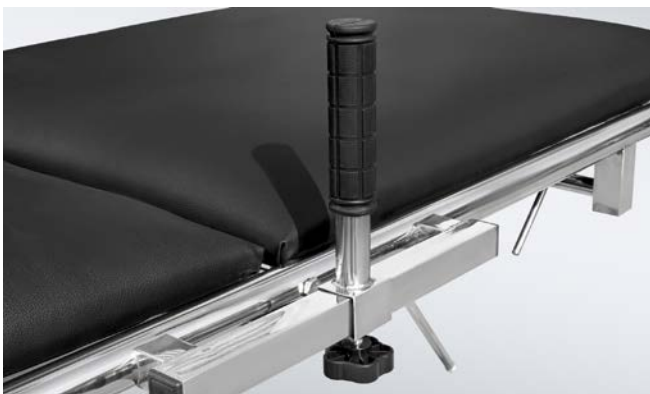

Handle

The handle can be moved on rails.



Leg support

Leg support can be adjusted in height and angle.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US