

# MC-D17

Electric Gynecological Table



*Better care  
Better comfort!*

# MC-D17

## ELECTRIC GYNECOLOGICAL TABLE

The humanized seat cushion heating design, multi-angle adjustment of leg support, and stable motor drive make the MC-D17 a beautiful and superior-performance gynecological examination table. The multi-angle adjustment and detachable functions of various parts make it meet various posture needs during gynecological examination.

### FEATURES


- Humanized design, seat heating function
- High-strength plastic cover, waterproof and antibacterial cushion, easy to clean and disinfect
- Electrically adjustable height and back section
- Leg section adjusted by gas spring
- The armrest automatically maintains a horizontal position,removable
- Remote control with lock button to prevent accidental touch
- The leg support can be adjusted at multiple angles to achieve the best production and inspection position
- 4 wheels for easy movement
- Removable full stainless steel filth basin

ITEM	PARAMETER	ITEM	PARAMETER
Overall size	1730*860*640~910mm	Lying surface size	1730*570mm
Back section up	0~85°±5°	Armrest size	560*100mm
Leg section down	0-80°±5°	Leg support swing out	360°
Height adjustment	640~910mm	Leg support height range	30cm
Back section length	830mm	Leg support forward and backward	45°,35°
Seat section length	470mm	Foot support height range	30cm
Leg section length	430mm	Safe working load	220kg

## PRODUCT DETAILS

### Foot support adjustment

Removable foot support horizontal adjustment range is 30cm.

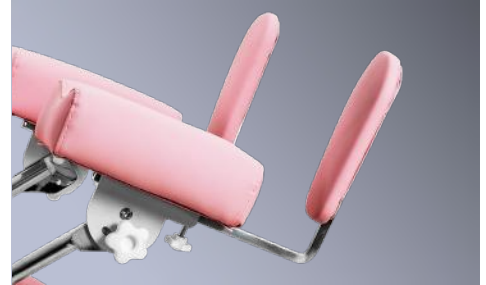

### Arm support adjustment

Removable arm support horizontal adjustment range is 30cm, can automatically maintain level during back adjustment.



### Filth basin

Removable full stainless steel waste basin.



### Pillow

Pillow can be adjusted and removed.



## PRODUCT DETAILS

### Hand remote

Bed adjustment by the hand remote,has lock button.




### Seat cushion

Removable seat cushion.



### Heat seating function

Heat seating function,humanized design.

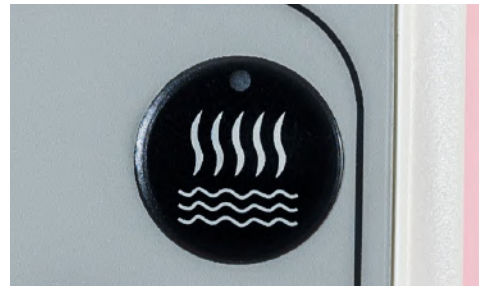

### Castors

Four dia 100mm castors with individual brake.



## ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**info@medik-medical.com | +86 139 6223 1292**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US