

# YA-EC-U03

## Electric Ultrasound Examination Table



*Better care  
Better comfort!*

YA-EC-U03

ELECTRIC ULTRASOUND EXAMINATION TABLE

DESCRIPTION

Automatic linen sheet change will improve the exam efficiency, especially for the ultrasonic room, and it will provide a good hygiene experience for patients, comes with very comprehensive position adjustment capabilities, and YA-EC-U03 which bears the features is an ideal option for an ultrasound exam room.

- MAIN FEATURES
- Steel with powder coating structure, durable, antibacterial and anti-corrosion
  - Seamless high-density mattress, antibacterial,corrosion-resistant and scratch-proof
  - Trendelenburg/Reverse Trendelenburg/height /Back and seat section adjustment function
  - Sheet roll holder on both ends of the bed
  - Automatic bed sheet change with automatic deviation correction function
  - Quick bed sheet changing system
  - Four castors ,dia 5 inches ,with central brakes for easy transport
  - Has emergency stop button
  - With battery backup for emergency situation

TECHNICAL  
PARAMETER

|                                     |                    |
|-------------------------------------|--------------------|
| Overall size                        | 2100*770*620-900mm |
| Height adjustment                   | 620-900mm          |
| Trendelenburg/Reverse Trendelenburg | 0~12°              |
| Back section                        | 0~70°±5°           |
| Seat section lifting                | 0~100mm            |
| Static working load                 | 250kg              |
| Dynamic working load                | 170kg              |
| Castor                              | Dia 5 inches       |
| N.W                                 | 145kg              |
| G.W                                 | 185kg              |

## PRODUCT DETAILS

### Pillow and Mattress

Ergonomically designed pillow for easy use,seamless high-density mattress,antibacterial,corrosion-resistant and scratch-proof ,easy to clean.



### Push handle

Push handle on the end of the bed ,easy to transport.



### Sheet roll holder

Sheet roll holder on both ends of the bed ,easy for automatically sheet change.



### Quick bed sheet changing device

There are quick bed sheet changing devices , making it easy to replace new sheets.



## Hand remote

The hand control is used to automatically change sheets, also for bed movements control.



---

## Height adjustment

Height adjusted by hand remote control from 620-900mm.



---

## Trendelenburg

Trendelenburg adjusted by hand remote control from 0~12 degree.



---

## Reverse Trendelenburg

Reverse trendelenburg adjusted by hand remote control from 0~12 degree.



---

## Sliding function

Bed surface can sliding backward and forward from 0~300mm by hand remote.



### Back section adjustment


Back section can be adjusted by hand remote from  $0\sim70^{\circ}\pm5^{\circ}$ .



---

### Seat section adjustment

Seat section can be adjusted by hand remote from 0~100mm.

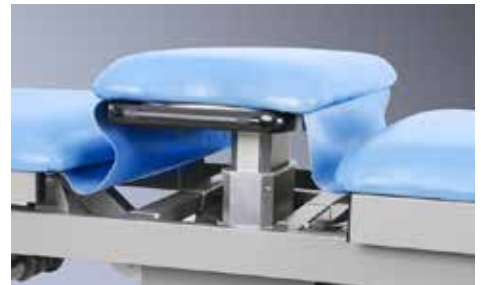

---

### Emergency stop button

Emergency stop button under foot section which can stop all movements urgently.



---

### Castors

Four pcs of dia 5 inches castors with central brakes at both sides of bed end.



## ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**info@medik-medical.com | +86 139 6223 1292**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US