

YA-EC-U06

Electric Gynecological Examination Table



*Better care
Better comfort!*

YA-EC-U06

/ELECTRIC GYNECOLOGICAL EXAMINATION TABLE

Automatic linen sheet change will improve the exam efficiency, especially for the ultrasonic room, and it will provide a good hygiene experience for the patient, the multi-angle adjustable leg rest can achieve a variety of inspection angles, all these features make YA-EC-U06 widely suitable for ultrasound rooms and gynecological examination rooms.

Main Features

- Steel with powder coating structure, durable, antibacterial and anti-corrosion
- Seamless high-density mattress, antibacterial, corrosion-resistant and scratch-proof
- The height and angle of leg rest can be adjusted
- Sheet roll holder on both ends of the bed
- Automatic bed sheet change with automatic deviation correction function
- Quick bed sheet changing system
- Four castors ,Dia 5 inch ,with central brakes
- Emergency stop button
- With battery backup for emergency situation
- With Filth basin

Technical parameter


Overall size	1680*1040*570-850mm
Height adjustment	570-850mm
Back section	0~65°±5°
Foot rest rotate angle	0~180°
Foot rest lifting angle	-20~20°
Filth basin adjusted distance	0~100mm
Static working load	250kg
Dynamic working load	170kg
Castor	Dia 5 inch
N.W	102kg
G.W	127kg

Main details



Height adjustment

Height adjusted by hand remote control from 570~850mm.

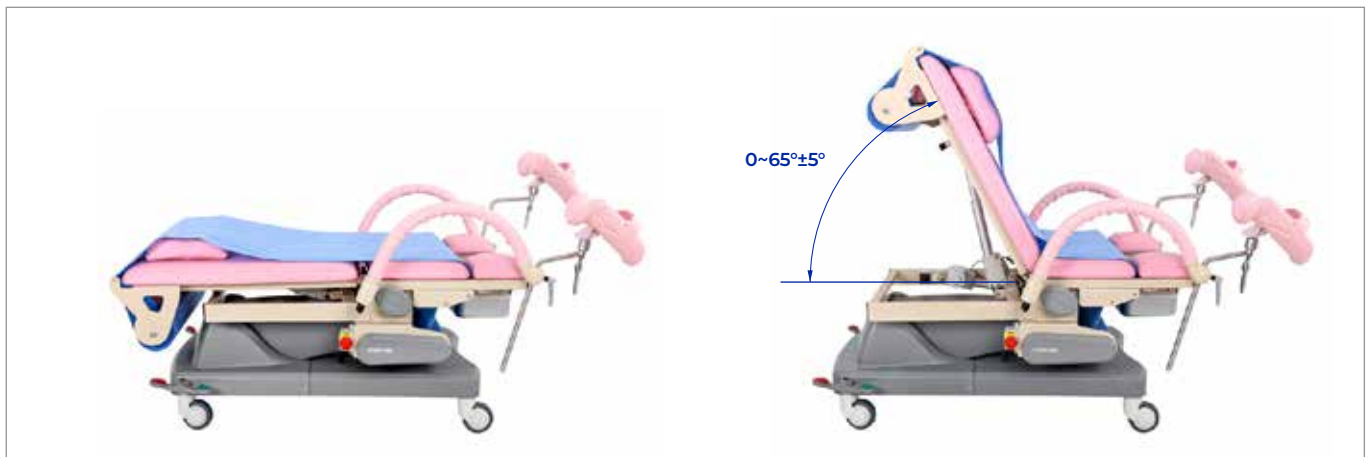

Back section adjustment

Back section can be adjusted by hand remote from 0~65°±5°.




Leg rest adjustment

The height and angle of leg rest can be adjusted through the locker and universal ball.



Pillow and Mattress

Ergonomically designed pillow for easy use, seamless high-density mattress, antibacterial, corrosion-resistant and scratch-proof, easy to clean.

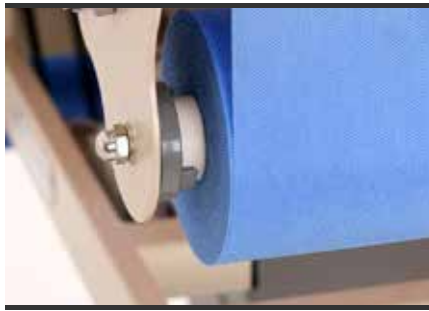

Sheet roll holder

Sheet roll holder on both ends of the bed, easy for automatically sheet change.



Quick bed sheet changing device

There are quick bed sheet changing devices, making it easy to replace new sheets.



Hand remote

The hand control is used to automatically change sheets, also for bed movements control.



Filth basin

Retractable filth basin under for easy use and clean.



Emergency stop button

Emergency stop button under foot section which can stop all movements urgently.



Castors

Four pcs of dia 5 inches castors with central brakes at both sides of bed end.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US