

MC-C03

Manual Gynecological Table



*Better care
Better comfort!*

MC-C03

MANUAL GYNECOLOGICAL TABLE

DESCRIPTION

The strong and durable stainless steel structure, thick waterproof and flame retardant cushions, and multiple position adjustments make this gynecological examination table not only meet the necessary examination needs, but also save costs.

MAIN FEATURES

- High grade stainless steel structure,durable and easy to clean
- Has trendelenburg /reverse trendelenburg position
- Waterproof and flame retardant cushions
- Leg section adjusted by gas spring
- Leg support can be adjusted at multiple angles to achieve the best production and inspection position
- Mounted on four units of anti-slip rubber feet
- Removable full stainless steel filth basin

TECHNICAL PARAMETER

Model	MC-C03
Lying surface size	1850*600*780mm
Back section up	0~60°±5°
Trendelenburg/Reverse-Trendelenburg	0~12°±2°
Mattress thickness	5cm
Safe working load	180kg
Package	Knock down design

PRODUCT DETAILS

Back section adjustment

Back section can be adjusted by gear from 0~60°.



Leg section adjustment

Leg section can be flipped under the seat section to achieve gynecological examination position.



Trendelenburg position

Trendelenburg position can be adjusted by latch from 0~12°.



Reverse trendelenburg position

Reverse trendelenburg position can be adjusted by latch from 0~12°.



PRODUCT DETAILS


Shoulder support

Shoulder support can be adjusted horizontally by knob.



Hand grip

Hand grip can be adjusted horizontally by knob.

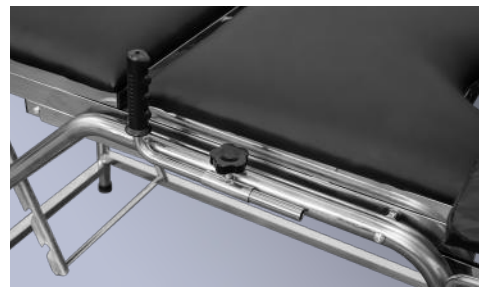

Leg support

Height of the leg support can be adjusted, the angle can be adjusted through the universal ball.



Filth basin


Removable full stainless steel waste basin.



PRODUCT DETAILS

Semicircular notch

Semicircular notch facilitates liquid flow out.

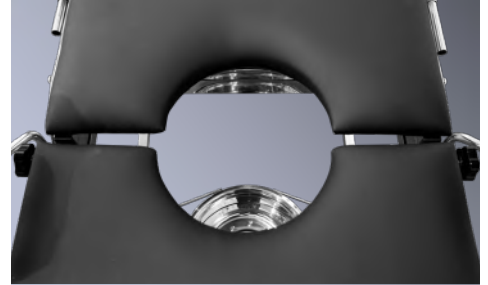

Rubber feet

Mounted on four units of rubber feet, anti-slip.



ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.