



DELIVERY BED
MODEL: MC-13

DELIVERY BED

Compact design, multifunctional design usage including delivery, gynecological examination, make this bed become a popular selection with a limited budget. Height adjustable is convenient for different height patient and medical staff, concealable leg section can convert the bed to a table position. All these features make MC-C13 as an ideal affordable product.

- Powder coating structure
- Cushioned with foam mattress
- Leather covered, waterproof
- Foot section mattress removable.
- Foot section mounted on silent rails sliding under the main frame.
- Adjustments of back rest by ratchet.
- Variable height adjustable by a hydraulic pump from 60 to 90 cm.
- Fixing clamp for anatomical knee crutches.
- Stainless steel sliding tray for liquid collection under the seat.
- Mounted on four rubber feet
- Knock down construction.

YA-MC-13 | Technical data

NO.	ITEM	PARAMETERS
01	Overall size (L*W*H)	180*80*60-90cm
02	Back section adjustment	0-70°
03	Loading capacity	170kg
04	Leg section adjustment	Pack up



Height adjust pedal

The height adjustment is controlled by the foot pedal, which is convenient for the doctor to directly control the height of the examination bed with his feet during the examination, which not only liberates the doctor's hands, but also saves effort.



Handle

The curved handle can adjust the front and back position, provide support during the delivery or examination of the patient, and can be stowed down when it is not needed.




Leg support

The leg support can be adjusted in height and curvature, which greatly improves the comfort of the patient.



Dirt Basin

Stainless steel dirt basin easy to clean, durable and replaceable.

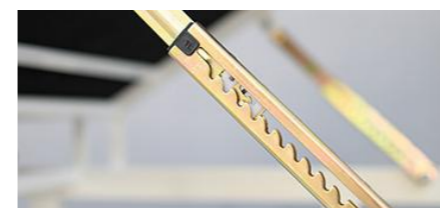

Back adjust

The adjustment of the back is controlled by saw-tooth, which can adjust each required angle. Pull the back to the highest point to level the back.



Bed legs

Rubber legs are installed on the bed legs to prevent noise from rubbing with the floor.



You Tube



Medik

Tel : +86 15250358182 / Fax : +86 51258688450

E-mail : info@medik-medical.com

Address : No.3 Jinxiu Road,Jinfeng Industrial Zone,Jinfeng
Town,Zhangjiagang,Jiangsu Province,China