

MC-D07

Delivery bed



*Better care
Better comfort!*

MC-D07

Delivery bed

DESCRIPTION

MC-D07 is a fully motorized and versatile LDR bed enables the obstetrician as well as patient through all stages of parturition, including labor, delivery & recovery stage with utmost comfort and ease. It can be easily adjusted to delivery position. Also electric controlled height and tilt adjustment greatly enhance the work efficiency of the nurses.

FEATURES

- Epoxy coated steel bed frame , with ABS engineering plastic side and head board
- The table surface is a waterproof PVC cushion with foam stuffed.
- Equipped with switch
- Ellescopic working table , to be controlled by linear bearing, telescopic length : 550mm
- Backrest can be adjusted electrically from 0 to 70°
- Height adjustable: 650mm-960mm
- Reverse-trendelenburg from 0 to 12°
- Trendelenburg 0 to 12°
- The mattress adopts high density memory sponge.
- With four 125mm castors with twin-wheel central braking system

STANDARD SPECIFICATION

Overall size	1950*950*650-960mm
Backrest adjustable	0-70°
Height adjustable	650-960mm
Reverse trendelenburg	12°
Trendelenburg	12°
Loading capacity	230kg



X-RAY TRANSLUCENT BEDBOARD

Full length X-ray translucent platform is equipped so X-ray photography can be done on the bed without moving the pregnant woman.



POSITION FIXER

The legs section and dirt basin can be fixed to prevent unnecessary slippage when adjusting the position.



CPR FUNCTION

Manual CPR handle on both the side of the back section. Release the back-rest for emergency situations.




DIRT BASIN

Large-capacity dirt basin, SUS 304 material is durable and not easy to rust.




REMOTE CONTROL

The hand control with clear icon is beneficial for medical staff which can control the bed position from a remote distance.

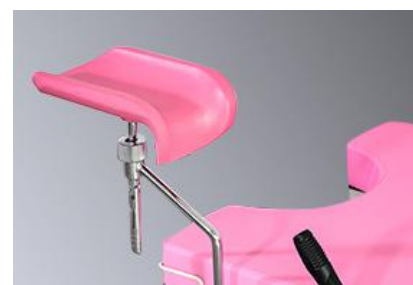

FOOT SUPPORT

The foot support is easy switch to be positioned to various positions. When it removed, the bed can be totally placed to be as a normal bed.

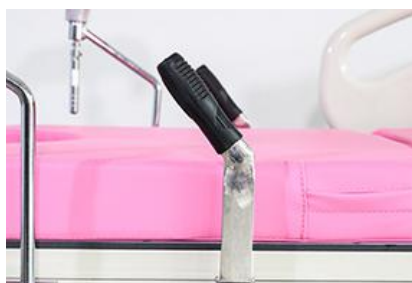

HANDLE

The curved handle can adjust the front and back position, provide support during the delivery or examination of the patient, and can be stowed down when it is not needed.



CENTRAL LOCK

Central lock system is convenient for the nurse to push the patient to the required position and quickly lock it with one key, which can appropriately reduce the workload of the medical staff.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US