

medik

Tel : +86 15250358182
E-mail : info@medik-medical.com

EMERGENCY CART

MK-C07

There is a trend to increase the color for the hospital interior area such as emergency or ICU room, so the patient will have a better feeling. The MK-C07 cart is an elegant colourful emergency trolley composed of 4 aluminum alloy and colored spray, equipped with accessories such as waste bin, oxygen cylinder, CPR board, etc. which can store medical supplies and facilitate medical staff to take them.

** Trolley designed to manage patient in emergency.*

** Structure: Sturdy, uni-body steel construction*

** Painted steel top rail / push handle*

** Removable plastic top*

** With I.V. transfusion hook at the top*

** With defibrillator shelf, cpr board, external power holder and plug board,*


oxygen tank holder, waste bin, sliding writing shelf, wire basket

** Four 100mm elegant castors and two with brakes*

** With 3 drawers and cupboard for storage.*

** Seal Lock or Centre Lock and Soft closing Drawer.*

SIZE : 680×480×1000MM ■
COLOR : WHITE/ RED/GREEN ■

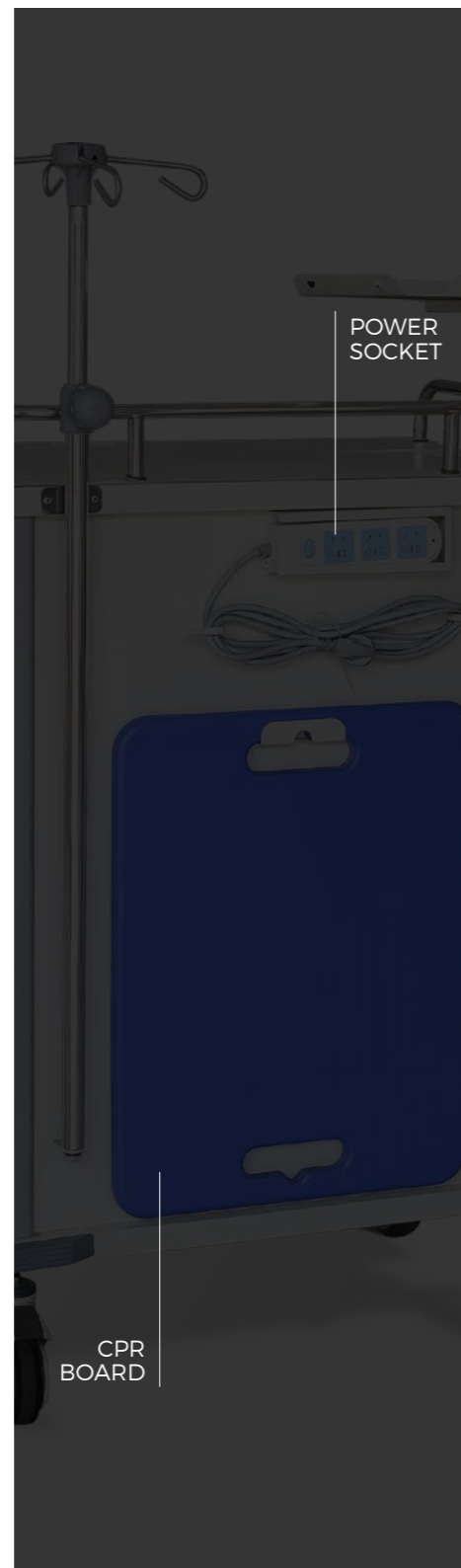

Extend work surface

on the left side of the cart and can be pulled out easily, which is convenient for medical staff to write



Top board

Stainless steel top board with stainless steel guardrails on three sides



Castors

Four 100mm elegant castors and two with brakes



Oxygen tank holder and waste bin

oxygen cylinder holder can hold 5L oxygen cylinders Waste bin on the left and right sides, which is convenient for medical staff to throw garbage



Drawer

3 same height drawers, each drawer has built-in divider to store items

