

YA-AS01

Ambulance stretcher



*Better care
Better comfort!*

YA-AS01

Ambulance stretcher

DESCRIPTION

The YA-AS01 ambulance stretcher with the removable stretcher surface offers excellent manoeuverability. Back section, leg section and height can be adjusted to meet different situation. It can flexibly transport the injured in various occasions such as gymnasiums, fields, and emergency rooms, light weight and convenient.

FEATURES

- These automatic loading ambulance stretchers are made from high-strength aluminum alloy.
- The foam cushion and the adjustable back, allow the maximum comfort.
- It is easy to operate, one person can operate and send the patient (with this stretcher) to the ambulance vehicle.
- The mechanical foldaway legs are controlled by the handles on one end of the stretcher.
- The legs can be folded under the stretcher automatically when loading onto the vehicle.
- Also the stretcher can be locked by the security devices on the ambulance vehicle.
- This type of stretcher can be characterized as light weight, durable, safe and easy to operate.


STANDARD SPECIFICATION

Material	High-strength aluminum alloy
Size(Low position)	1960*550*340mm
Size(High position)	1960*550*980mm
Backrest adjustable	0-65°
Footrest adjustable	0-40°
Load capacity	160kg



Foot section adjust

The foot section is adjusted by the gas spring . It can be adjusted according to the condition of the injured person.

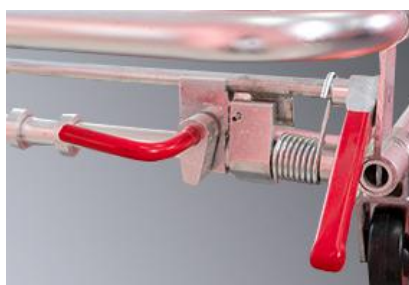

Safety lock

On the side of the Legs retracted handle, there is a safety lock. Before retracting the legs, the lock needs to be opened. This mechanism prevents accidental touch.



Stretcher surface separation

By pulling the handle, the bottom stretcher surface and the base can be separated to facilitate patient transfer.




Legs retracted

When pushing the stretcher to the ambulance, you can press the handle of the head to help retract the legs, which is quick and convenient to transfer the patient to the ambulance. At the same time, the height position can also be adjusted by pulling the handle.



Height adjust

Pull the handle at the same time on the head and legs of the stretcher to adjust to the lowest position.

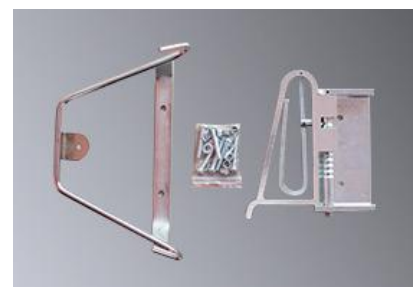

Fixing device

Each ambulance transporter comes with a fixing device that can be attached to the ambulance to prevent slipping.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US