

YA-AS08

Ambulance stretcher



*Better care
Better comfort!*

YA-AS08

Ambulance stretcher

DESCRIPTION


The YA-AS08 ambulance stretcher with the removable stretcher surface offers excellent manoeuvrability. Back section, leg section and height can be adjusted to meet different needs. It can flexibly transport the injured in various occasions such as gymnasiums, fields, and emergency rooms.

FEATURES

- The surface is made of foamed cushion, vinyl protective cover that is durable and easy to clean
- Infusion stand optional
- Folding legs controlled by handles on both sides of the stretcher
- The stretcher can be locked by the fixing devices when on the ambulance car
- Trendelenburg
- Reverse Trendelenburg
- The height is adjustable


STANDARD SPECIFICATION

Material	High-strength aluminum alloy
Size(Low position)	1960*550*340mm
Size(High position)	1960*550*980mm
Backrest adjustable	0-65°
Footrest adjustable	0-40°
Trendelenburg	0-30°
Reverse Trendelenburg	0-30°
Load capacity	160kg

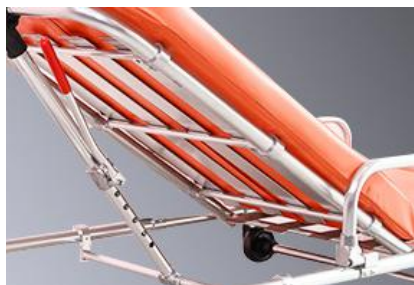

FOOT SECTION ADJUST

The foot section is adjusted by the gas spring . It can be adjusted according to the condition of the injured person.

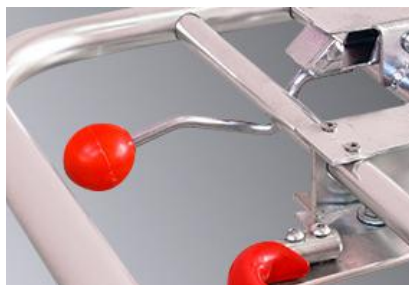

STRETCHER SURFACE SEPARATION

By pulling the handle, the bottom stretcher surface and the base can be separated to facilitate patient transfer.



LEGS RETRACTED

When pushing the stretcher to the ambulance, you can press the handle of the head to help retract the legs, which is quick and convenient to transfer the patient to the ambulance. At the same time, the height and Trendelenburg & Reverse Trendelenburg position can also be adjusted by pulling the handle.



HEIGHT ADJUST

Pull the handle at the same time on the head and legs of the stretcher to adjust to six different heights.



TRENDELEBURG POSITION

By pulling the handle of the head, the head can be brought down quickly, reaching trendelenburg position.



REVERSE TRENDELEBURG POSITION

By pulling the handle of the leg section, the legs can be lowered at five different angles to suit any emergency situation.



AUXILIARY WHEEL

There are spare wheels on the six positions of the stretcher, which can be moved by the wheels even after the stretcher is pushed onto the ambulance, which saves effort.



FIXING DEVICE

Each ambulance transporter comes with a fixing device that can be attached to the ambulance to prevent slipping.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US