

# YA-D5-13

New ICU bed



*Better care  
Better comfort!*

# YA-D5-13

## New ICU bed

### DESCRIPTION

The YA-D5-13 bed's numerous features and ergonomic design provide safer patient mobility and convenience for caregiver's operation. Compact and sturdy structure, 5 functions with manual CPR, which can meet the needs of ICU patient care.

### FEATURES

- With embedded railing control, visible icons, with nurse control panel facing outside and limited button facing the patient.
- Central brake system is adopted with dual braking pedals on both sides for easy manipulation
- Bumper caster is located on four corners of the bed
- Manual CPR lever is located on both sides
- Angle indicator is achieved on the railing, especially the 30 degrees for back section adjustment is clearly marked for frequent use.

### STANDARD SPECIFICATION

Model	YA-D5-13
Overall size	2230*1000*465~765mm
Bed surface	1937*865mm
Backrest adjustable	0-70°
Footrest adjustable	0-25°
Height adjustable	465-765mm
Trendelenburg	0-12°
Reverse trendelenburg	0-12°
Loading capacity	230kg



**EXTERNAL SIDE RAIL CONTROL**

The clear icon of the hand control is beneficial to the operation of staff



**CASTERS WITH CENTRAL BRAKE**

Central brake located at the foot of the bed (Red for locking, and another for unlocking)



**BUMPER WITH IV POLE HOLDER**

Bumper around the bed, with IV pole holes on the bumper.



**FOOTREST**

Feet section can be adjusted by racks, which can raise level the feet, promote blood circulation.



**INTERNAL SIDE RAIL CONTROL**

Internal side rail control patient can only control back and leg adjustment.



**CPR FUNCTION**

Release the back-rest for emergency CPR.




**ANGLE INDICATOR**

Easy-to-read angle indicators located on the head side rails highlight elevation in the 30°-45° range, assisting caregivers in.



**POWER SWITCH**

Power switch located at the head side to avoid the mis-operation.

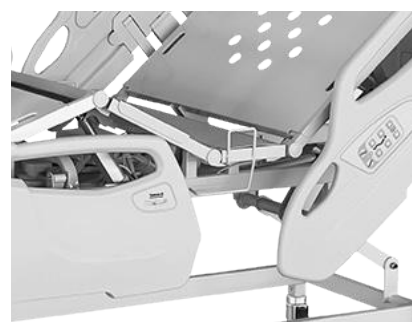

**EASY TO GRIP SIDE RAILS**

Ergonomically designed side rails to protect patients and can easily grip them for egress the bed.

## ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**info@medik-medical.com | +86 199 6269 8550**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US