

# YA-DH5-2

Home Care Bed



*Better care  
Better comfort!*

# YA-DH5-2

# Home Care Bed

## DESCRIPTION

YA-DH5-2 is an economical electric home care bed with multi-functions that can be adjusted to meet the different needs of peoples. The full-length side rail can prevent peoples from falling risk. This bed significantly reduces costs and meets the low budget needs of home care and nursing care center.

## FEATURES

- Electrically operated profiling bed for community and nursing home use
- Safe work load (SWL) 180Kg
- Variable height (350mm -750mm)
- Adjustable self-locking side rails providing additional safety which can be released at both head and foot ends
- Remote handset control
- High-quality castors (castor diameter: 100mm)
- Quick and easy assembly in limited space

## STANDARD SPECIFICATION

Model	YA-DH5-2
Overall size	2140*1010 *350-750mm
Lying size	1900*900 mm
Backrest adjustable	0-75°
Footrest adjustable	0-40°
Height adjustable	350-750mm
Trendelenburg	0-12°
Reverse-Trendelenburg	0-12°
Loading capacity	180kg



**CONTROL SYSTEM**

Total 4pcs linear heavy duty type actuator will guarantee different position adjustment and safety loading 180kg.



**BED PLATFORM**

Sturdy metal welded platform surface with mattress holder, 4 section board.




**SIDE RAILING STOPPER**

There are stoppers under the headboard and footboard, which can prevent the side rail from sliding down and prevent safety accidents.



**HAND PADDLE**

The height of the full-length guardrail can be adjusted by the wave hand paddle.

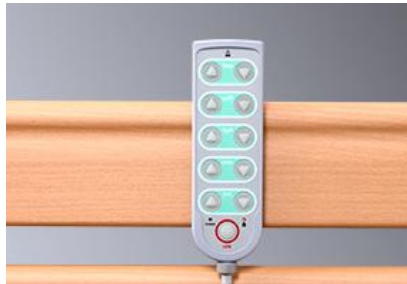

**REMOTE CONTROL**

The hand control with clear icon is beneficial for medical staff which can control the bed position from a remote distance.



**CALF BOARD**

Calf board can be adjusted by racks, which can raise level the feet, promote blood circulation if necessary.



**CASTOR**

Four  $\phi$ 10 cm swivel casters with brakes, easy to move the patient and fixed in any desired position.



**HEAD AND FOOT BOARD**

Simple head and foot boards, the top cover is plastic, and the 2 motors are exposed outside, which greatly reduces the cost.

## ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**info@medik-medical.com | +86 139 6223 1292**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US