

YA-DH8-1

Winglet lying surface hospital
patient bed



*Better care
Better comfort!*

YA-DH8-1

WINGLET LYING SURFACE HOSPITAL PATIENT BED

DESCRIPTION

YA-DH8-1 is a medical bed designed with future concept, aluminum casting frame makes the bed for long -term use and easy for disinfection. The unique butterfly sleeping platform provides the patient special comfort which can not be experienced from any kind of a single mattress.

FEATURES

- Whole aluminum alloy frame
- Unique telescopic aluminum railing foldable and can be sliding under to the platform
- Unique butterfly sleeping platform can distribute the pressure of the body evenly
- Manual CPR lever located on both sides of the bed
- Bumper caster located on four corners
- Central brake system with 150mm diameter
- Auto-counter to provide the patient from sliding onto foot end
- Hand remote control can stuck on the railing
- Detachable head/foot board with easy lock
- Bed extension in mechanical way (optional)
- Auto-regress to remove the pressure of the abdomen when adjust the backrest (optional)
- Retractable 5th wheel for direction(optional)

STANDARD SPECIFICATION

Overall size	2103*990*465-790mm
Back section up/down	0-70°
Foot section up/down	0-35°
Height adjustment	465-790mm
Trendelenburg	0-12°
Revere- Trendelenburg	0-12°
Loading capacity	240kg
Castor	Dia 150mm



Lying surface

Unique butterfly platform can distribute the patient weight evenly from each single butterfly knot. Thus it can provide great sleeping experience and can minimize the risk of bedsores.




Section adjustment

Height, back and foot section can be adjusted by remote controller, auto-counter to provide the patient from sliding into foot end.



Telescopic side rail

Unique telescopic aluminum railing foldable and can be sliding under to the platform, protect patient safety and make it easier for patients to get out of bed.

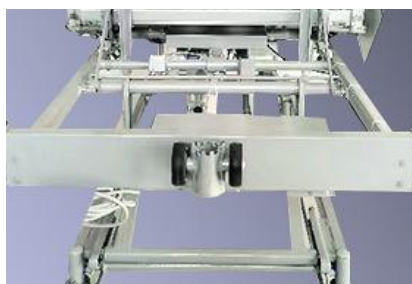

Removable head and foot board

The headboard and footboard of the bed are removable, increasing the flexibility of bed use.



Bumper caster

Bumper caster located on four corners of the bed to prevent the bed from being damaged.



Emergency stop button

Emergency stop buttons are located on both sides of the bed and can emergency stop the bed's progress.



CPR function

Manual CPR lever on both sides of the bed, also hand remote control has electric CPR function.



Castors

4 luxurious silent castors with central brake system.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US