

medik



YA-DS-M03

& MANUAL
DIALYSIS CHAIR

**BETTER CARE
BETTER COMFORT**

Zhangjiagang
Medi medical Equipment
Co., Ltd



Tel : +86 15250358182 / Fax : +86 51258688450
E-mail : info@medik-medical.com
Address : No.3 Jinxu Road,Jinfeng Industrial Zone,Jinfeng
Town,Zhangjiagang,Jiangsu Province,China

YA-DS-M03 MANUAL DIALYSIS CHAIR

YA-DS-M03 is a medical chair used widely in different hospital areas, such as blood donor, dialysis or infusion thanks to its functional design. The Ergonomic design provides both patients and nurses with comfortable experience.

- * Chair designed according to ergonomic human engineering
- * With comfortable mattress, PU cover, waterproof
- * Adjustable by gas spring
- * Four casters and two with brake



YA-DS-M03 DETAILS DISPLAY



SEAT ADJUST

The seat section can be adjust by the handle on the left side ,with different angles to provide more comfortable experience to the users.



BACK/FOOT SECTION ADJUST

The handles on the right side control the angle of back and foot section . Toggle the handle to control the back and legs and adjust it to a comfortable angle for the user.




IV POLE HOLDER

High-quality epoxy-painted material is used to place the infusion stand, and there are knobs to fix it.



PILLOW

High-density foam pillow can be easily attached and removable to the head section for comfort.

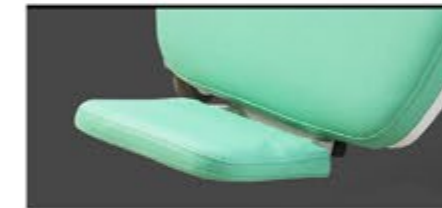

FOOTREST

Foldable footrest can give foot support during a semi fowler position and prevent the patient from sliding downwards.



CASTOR

5" casters with individual brakes allow for better handling and control.



Technical parameter

01	Overall size	1880*830*560mm
02	Back adjust	-10° to 80°
03	Seat adjust	0° to 30°
04	Leg adjust	-80° to 10°
05	Load capacity	200KG
06	Back size	190*490*120mm