

BETTER CARE
BETTER COMFORT

Zhangjiagang
Medi medical Equipment
Co., Ltd



YA-DS-M05 MANUAL DIALYSIS CHAIR

YA-DS-M05 is a medical chair used widely in different hospital areas, such as blood donor, dialysis or infusion thanks to its functional design. The Ergonomic design provides both patients and nurses with comfortable experience.

- * Ergonomic design , makes the whole process easier and enjoyable.
- * Back and leg section adjustable by manual switch.
- * Powder coated steel chair frame
- * High-quality medical polyurethane foam molding
- * Anti-static and anti-fouling properties



YA-DS-M05 DETAILS DISPLAY



BACK/FOOT SECTION ADJUST

The handles on the right side control the angle of back and foot section , Toggle the handle to control the back and legs and adjust it to a comfortable angle for the user.



ARM SWING OUT

armrest can swing out 0-180°




LEATHER AND FRAME

comfortable genuine leather and cover and high quality powder coated



IV POLE HOLDER

High-quality epoxy-painted material is used to place the infusion stand, and there are knobs to fix it.

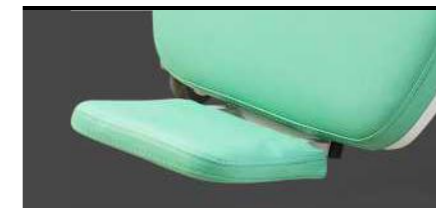

FOOTREST OPTIONAL

Foldable footrest can give foot support and prevent the patient from sliding downwards.



CASTOR OPTIONAL

5" casters with individual brakes allow for better handling and control.



Technical parameter

01	Overall size	1900*500*550mm
02	Back adjust	0° to 90°
03	Leg adjust	0° to 90°
04	Armrest swing out	0° to 180°
05	Armrest height adjustable	100mm
06	Load capacity	200kg

Accessory: IV Pole;
Paper roll holder;
Oxygen cylinder holder;
Dining board