



Tel : +86 15250358182
E-mail : info@medik-medical.com

RECLINING PHLEBOTOMY CHAIR

YA-DS-R03

YA-DS-R03 Chair Recliner offers comfort and durability to patients in hospitals, clinics or assisted living facilities. This chair has three traditional positions - a deep recline, elevated leg rest and full upright, operated by a retractable lock bar, securing the chair in each desired position.

- Ideal for long term care facilities, retirement homes and nursing homes.
- Three positions: upright, deep recline, and elevated leg rest/footrest
- Retractable lock bar secures the chair in the desired position
- Padded push-handle for easy maneuverability.
- Extra-wide padded armrests and side panels maximize comfort
- Blow-molded tray locks in 5 positions and can be stored on the side of the chair
- Durable epoxy coated welded steel frame is easy to clean and maintain
- Articulated headrest quickly adjusts to maximize comfort
- Chair is easy for patients to move and reposition
- Four swivel casters with locking casters in rear



ITEM	PARAMETERS
Model	YA-DS-R03
Overall size(H*W*D)	1200 x 635 x 1040mm
Seat size(H*W*D)	530 x 482x 482mm
Tray size	650*390mm
Floor to Seat Height	660mm
Width between arms	508mm
Loading capacity	113kg




Retractable lock bar

Retractable lock bar has 3 different adjustments with 3 different positions (Upright position, deep recline position, elevated leg rest position) can meet the different needs of patients



Footplate

Allow patients to put their feet onto it for transferring the patient

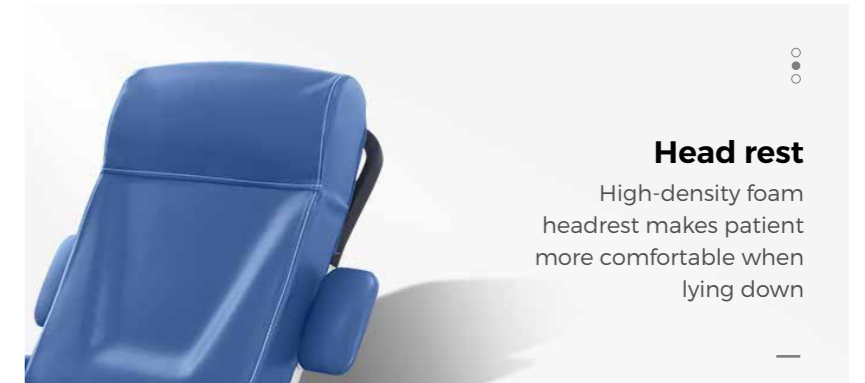

Head rest

High-density foam headrest makes patient more comfortable when lying down



Side panel

Large, blow-molded tray locks in 5 positions for different size patient and can be stored on the side of the chair



Casters

5"Total casters with individual lock allow for better handling and control.