

YA-SH3-1

Semi Electrical Home
Care Hospital Bed



*Better care
Better comfort!*

YA-SH3-1

SEMI ELECTRICAL
HOME CARE HOSPITAL BED

DESCRIPTION

YA-SH3-1 is the ideal bed for home care,we designed the most often used functions to be electric, and the less often used functions to be manual,this will maximizes care efficiency and minimize the cost for users.It’s durable quality and appearance close to home decoration can bring the best home care experience.

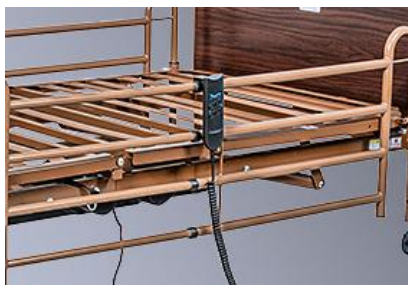
- FEATURES
- Made of steel frame with coating surface
 - Back and foot control by motor, Hi/Lo control by hand crank
 - The plastic grain panels are break and scratch resistant
 - Channel frame construction provides superior strength and reduced weight.
 - Zinc-coated spring deck
 - Unique motor is completely self-contained to reduce weight and noise
 - The motor assembly can be installed and removed with patient in bed
 - Bed ships in two cartons
 - Four 75mm diameter castors and two of them with brake.

STANDARD
SPECIFICATION

Overall size	2185 x 905mm
Mattress size	2015 x 880mm
Adjustable height	385-585mm
Angle of back section	0 – 70°
Angle of leg section	0 – 35°
Net Weight	72.4kg
castors	Diameter 75mm


Section adjustment

Back section and foot section can be adjusted by remote hand control separately, also there is the auto-counter function which can adjust patient back board and leg board together easily.


Side rail

The height of side rail can be adjusted manually and can be auto-retract, it can ensure patient safety during the position adjustment process.


Height adjustment

Height of the bed can be adjusted by manual crank on bed end.


Remote hand control

Easily identifiable icons on the remote, and when you press the two half circles on the middle, the back and foot section can be raised or lowered at the same time.


Bed frame

The steel channel structure can make the bed strong and light, so it is very convenient to pack and ship, also equipped with the mattress fixtures.


Motors

Equipped with well-known brand boat type motor, noiseless but high performance.


Castors

Four pcs dia 75mm diameter noiseless double side castors on the bottom and two of them with brake.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US